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(Válida únicamente en los EE.UU.)

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IMPORTANT SAFEGUARDS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING
DANGER- TO REDUCE THE RISK OF ELECTROCUTION:

- ALWAYS unplug this product immediately after using.
- DO NOT handle the plug with wet hands.
- DO NOT place or store the unit where it can fall or be pulled into a tub, sink or toilet.
- DO NOT reach for a unit that has fallen into water. Unplug immediately.
- DO NOT use while bathing.
- DO NOT place in or drop into water or any other liquid.

WARNING - TO REDUCE THE RISK OF BURNS, ELECTROCUTION, FIRE OR INJURY TO PERSONS:

- A product should never be left unattended when plugged in.
- DO NOT excessively bend or pull the cord. Unplug the cord by grasping the plug—not the cord.
- Keep the cord away from heated surfaces.
- Use this product only for its intended use as described in this manual. DO NOT use the product for any other purpose.
- DO NOT use attachments not recommended by HoMedics.

PowerDent Oral Irrigator is perfect for cleaning around braces, crowns or bridgework.

This unit can be used with water or your favorite brand of mouthwash.

- DO NOT use outdoors or operate where aerosol products are being used or where oxygen is being administered.
- NEVER insert irrigator-tip into or remove it from the universal handle while the product is on.
- NEVER operate the product if it has a damaged cord or plug, if it is not working properly, if the charging base has been dropped or damaged, or if it has been immersed in water. Return the product to HoMedics service center for examination and repair.
- This product should never be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- This unit should not be used by children or invalids without close supervision.
- NEVER block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- NEVER use while sleeping or drowsy.
- NEVER drop or insert any object into any opening or hose.
- Not recommended for children under 3.
- For household use only.
- Fill reservoir with water only unless otherwise specified by manufacturer.

SAVE THESE INSTRUCTIONS

PowerDent™ Oral Irrigator

⑧ 4 COLOR CODED IRRIGATOR TIPS

HOW TO USE POWERDENT ORAL IRRIGATOR


This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not change the plug in anyway.

1. Plug the cord into any standard electrical outlet. If a wall switch controls the outlet, be sure it is turned on.
2. Remove the reservoir. When removing reservoir use one hand to firmly hold the base ① of the unit in place, and with the other hand, lift the reservoir straight up. Fill with lukewarm water. Reinstall the reservoir onto the unit.

Note: Prior to initial use, rinse the unit by letting a full reservoir of water pass through it. If the unit is used with a mouthwash solution, it should be rinsed with clean water after each use to remove any traces of mouthwash. If the unit has not been used for an extended period of time, to maintain hygiene it is recommended that you follow the directions for initial use, letting a full reservoir of water pass through.

3. Select an irrigator tip ⑧ for your personal use (tips are color-coded for multiple users). Lift the irrigator handle ③ out of the base and attach the irrigator tip to the handle so that the colored tip is facing the same direction as the water stop switch (Figure 1).
4. Adjust the pressure control on the base unit to minimum pressure (turn clockwise) for first-time use. Gradually increase the pressure over time until you find the setting most comfortable to you.

Figure 1

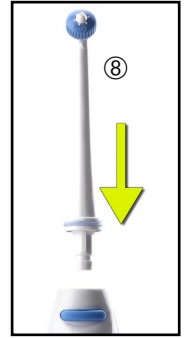Note: PowerDent Oral Irrigator comes with 4 color-coded irrigator tips for multiple users.

Figure 2

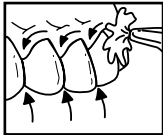
Figure 3

- Before turning the unit on, check that the water stop switch ⑤ on the handle is in the off position (away from the tip). After turning the unit on by pressing the on/off button ⑩ on the base unit, lean over the sink and position the tip in your mouth. Move the water stop switch to on (towards the tip) to begin the cleaning and massage process.

Note: Keep your mouth slightly open and hold your elbow at mouth level so that the water can run back into the sink.

- Direct the flow along the gumline (Figure 2) and also between your teeth. Do not point the irrigator tip directly at your gums. For best results, begin with your back teeth and work forward. Consult your dental professional for further recommendations on use.
- To adjust the flow of the water from a stream to the hydromassage spray (or spray to stream), rotate the colored irrigator tip (Figure 3). Use the stream for more pinpoint cleaning and the hydromassage spray for a more general cleaning and massage.

Note: To avoid splashing, adjust the flow from stream to spray while the unit is out of your mouth, over the sink and pointed down into the sink.

- When finished, move the water stop switch to the off position and turn off the unit. Remove the irrigator tip by pressing the release button ⑥ and replace in the covered storage. Empty any liquid remaining in the unit. Unplug the unit.

Note: NEVER leave the unit "On" when not in use. Water pressure will continue to build and may cause damage to the unit. You must turn the water stop switch ⑤ and the On/Off Button ⑩ to the Off position when not in use.

Note: We recommend soaking the irrigator tip in a denture cleaner (commercially available) every 2 to 4 weeks to keep it hygienically clean.

Caution: All servicing of this unit must be performed by authorized HoMedics service personnel only.

CLEANING AND MAINTENANCE

- After use, empty any liquid remaining in the unit.
Note: If the unit is used with a mouthwash solution, a full reservoir of water should be run through after each use to rinse it and remove any traces of mouthwash.
- Rinse the irrigator tip and the handle separately for several seconds under running water. Note: Be sure to point the handle downward (upside down) when rinsing. Do not allow water to drain into the irrigator handle.
- After cleaning, place the irrigator tip back in the covered storage and the handle on the base.
- Unplug the unit.
- Clean the base periodically with a soft cloth or paper towel. The base must be dry before plugging it into an electrical outlet.

MEDICAL SAFEGUARDS

Please note the following safeguards for your continued safety.

- Similar to starting any new oral hygiene product, contact your dentist before using PowerDent Oral Irrigator if you had oral or periodontal surgery.
- Contact your dentist if excessive bleeding occurs or if bleeding continues after 2 weeks.

Caution: Opening the handle will destroy the appliance and invalidate the warranty.

Note: During the first few days of using, your gums may bleed slightly. Any bleeding should stop within a few days. If your gums continue to bleed for a period of 2 weeks or more, discontinue use and consult your dentist.

ONE YEAR LIMITED WARRANTY

(Valid in USA only)

HoMedics, Inc., guarantees this product free from defects in material and workmanship for a period of one year from the date of purchase, except as noted below.

This HoMedics product warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of HoMedics. This warranty is effective only if the product is purchased and operated in the USA. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under warranty. HoMedics shall not be responsible for any type of incidental, consequential or special damages. All implied warranties, including but not limited to those implied warranties of fitness and merchantability, are limited in the total duration of one year from the original purchase date.

To obtain warranty service on your HoMedics product, either hand deliver or mail the unit and your dated sales receipt (as proof of purchase), postpaid, along with check or money order in the amount of \$5.00 payable to HoMedics, Inc. to cover handling.

Upon receipt, HoMedics will repair or replace, as appropriate, your product and return it to you, postpaid. Warranty is solely through HoMedics Service Center. Service of this product by anyone other than HoMedics Service Center voids warranty.

This warranty provides you with specific legal rights. You may have additional rights which may vary from state to state. Because of individual state regulations, some of the above limitations and exclusions may not apply to you.

For more information regarding our product line in the USA, please visit: www.homedics.com or www.turnonyoursmile.com

Due to continuous product improvements the pictures may differ slightly from the actual product.

Mail to:
HoMedics
Service Center Dept. 168
3000 Pontiac Trail
Commerce Township, MI
48390

e-mail:
cservice@homedics.com

© 2003 HoMedics, Inc. and its affiliated companies. All rights reserved. HoMedics® is a registered trademark of HoMedics, Inc. and its affiliated companies. PowerDent™ is a trademark of HoMedics, Inc. and its affiliated companies. All rights reserved.

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