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(Válida sólo en los EE.UU.)

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IB-HD270C

PowerDent™ ultra with Smart Sensor™ Technology

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IMPORTANT SAFEGUARDS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

Note: When the rechargeable cells are running empty, the toothbrush automatically turns off. This protects the rechargeable cells from deep discharging and signals that the power unit should be charged.

READ ALL INSTRUCTIONS BEFORE USING DANGER —TO REDUCE THE RISK OF ELECTROCUTION:

- DO NOT handle the plug with wet hands.
- DO NOT place or store the charger base where it can fall or be pulled into a tub, sink or toilet.
- DO NOT reach for a charger base that has fallen into water. Unplug immediately.
- DO NOT use while bathing.
- DO NOT submerge the charger base in water or any other liquid.
- This charger base is designed for use in 120-volt, 60 Hz outlet. NEVER force the plug into an outlet.

WARNING—TO REDUCE THE RISK OF BURNS, ELECTROCUTION, FIRE OR INJURY:

- DO NOT excessively bend or pull the cord. Unplug the cord by grasping the plug—not the cord.
- Keep the cord away from heated surfaces.
- Use this product only for its intended use as described in this manual. DO NOT use the product for any other purpose.

- DO NOT use attachments not recommended by HoMedics.
- DO NOT use outdoors or operate where aerosol products are being used or where oxygen is being administered.
- NEVER insert brush head into or remove it from the universal handle while the product is on.
- NEVER operate the product if it has a damaged cord or plug, if it is not working properly, if the charging base has been dropped or damaged, or if it has been immersed in water. Return the product to HoMedics service center for examination and repair.
- This product should never be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- This unit should not be used by children without adult supervision.
- Not recommended for children under 3.
- Avoid using a brush with crushed or bent bristles.
- NEVER use while sleeping or when drowsy.
- For household use only.

SAVE THESE INSTRUCTIONS

Note: The handle may warm up during charging. This has no effect on the unit's operation.

PowerDent Ultra comes with 3 color coded brush heads for multiple users.

Caution: All servicing of this unit must be performed by authorized HoMedics service personnel only.

Dentists recommend that you brush your teeth gently to avoid tooth and gum abrasion. Exclusive HoMedics' Smart Sensor™ technology activates the brushing action when just the lightest amount of pressure is applied, to help protect teeth and gums and reinforce proper brushing habits.

CHARGING INSTRUCTIONS

Follow this procedure before you use the PowerDent Ultra for the first time.

To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

1. Plug the charger base into a 120-volt outlet. Place the handle into the base and allow it to charge for at least 16 hours before using it for the first time.
2. The charger indicator light on the handle will glow red when the tooth brush is correctly placed in the charger base. This is an indication that the handle is charging. The indicator will turn off when the toothbrush is removed from the base.
3. The optimum voltage for charging is 120V.
4. To maintain the unit's full performance, it should always be replaced back on the charging base after use.

Note: Charge for at least 16 hours before using the unit for the first time.


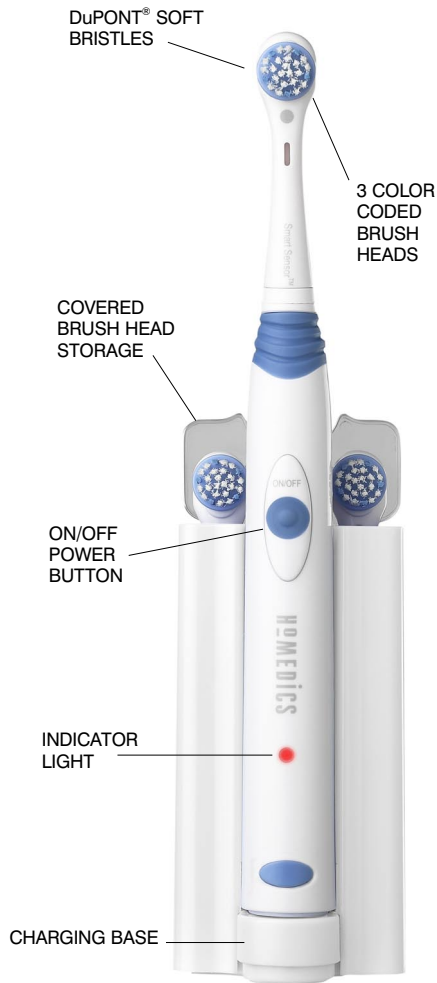
Figure A

Figure B

Figure C

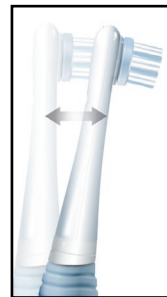HOW TO USE POWERDENT ULTRA

1. Attach the brush head to the handle by positioning the brush head at a 90-degree angle to the left of the handle (See Figure A). Turn the brush head counter clockwise one-quarter turn, so the bristles are facing the same direction as the ON/OFF button, and the brush head will snap into place. The Smart Sensor logo will be facing forward as will the bristles. (See Figure B)
2. Before brushing, wet the brush head to moisten, then apply toothpaste.
3. **TO USE WITH SMART SENSOR TECHNOLOGY**


Push the ON/OFF button briefly, and the corresponding LED light will glow green, indicating that the power is ON. Apply brush head against your teeth to activate the toothbrush (See Figure C). **Unit will not operate (no function) until pressure is applied to the brush head.** If you do not apply pressure to the brush head for 15 seconds the unit will turn off, along with the green LED light.

TO USE WITHOUT SMART SENSOR TECHNOLOGY


Press and hold the ON/OFF button for about 2 seconds, after which time the unit will automatically begin working. This option, to override the Smart Sensor feature, is available for those who have learned how to brush with a light touch or would prefer not to use this feature.

4. Guide the brush head to your teeth, placing it at a 45 degree angle towards the gum line. Gently guide the brush head from tooth to tooth, spending a few seconds on the surface of each tooth.


Outside teeth

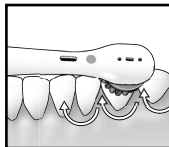Inside teeth

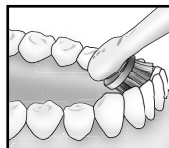Chewing surface

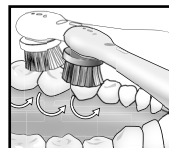Note: We recommend soaking the brush head in a denture cleaner (commercially available) every 2 to 4 weeks to keep it hygienically clean.

- Brush the gum lines as well as the teeth, starting with the outside teeth and gums, followed by the insides and finish with the chewing surfaces. Ensure that you spend equal time in all 4 quadrants of your mouth.
- Do not press too hard or scrub with this unit. Remember the unit is activated with the lightest pressure to help protect teeth and gums.
- The auto timer feature will automatically turn the unit off after the 2-minute dentist recommended brushing time. Note: To turn the unit off before the 2 minutes have elapsed, simply press the ON/OFF button. Ensure the unit is off when traveling with the PowerDent Ultra.
- Rinse mouth thoroughly with water.
- Remove brush head and rinse under running water.
- Brush heads should be replaced every 2-3 months to ensure optimal performance.

HoMedics replacement brush heads (RH-P) for this model are available at your local retailer or check www.homedics.com for a list of retailers that offer these replacement heads.

CLEANING AND MAINTENANCE

- After use, turn the toothbrush off and rinse the brush head for several seconds under running water.
- Take the brush head off the handle and rinse both parts under running water and wipe them dry. Note: Be sure to point the handle downward (upside down) when rinsing. Do not allow water to drain down into handle.
- Place the handle back on the charging base. Note: For optimal charging base stability, run the power cord through the right or left slot on the back of the charging base.

- When the unit is not in use for an extended period, unplug the cord and store.
- Clean the charger base periodically with a soft cloth or paper towel. The charger base must be dry before plugging it into an electrical outlet.

MEDICAL SAFEGUARDS

Please note the following instructions for your continued safety.

- Similar to starting any new oral hygiene regime, contact your dentist before using PowerDent Ultra if you have had oral or periodontal surgery.
- Contact your dentist if excessive bleeding occurs or if bleeding continues after 2 weeks.
- This product is designed for cleaning teeth and gums only. Do not use for any other purpose.

SPECIAL INSTRUCTIONS

REMOVING AND DISPOSING OF NI-CAD BATTERIES

This unit is equipped with long-lasting, rechargeable nickel-cadmium batteries. To help protect our environment, they should be removed at the end of the unit's life for safe disposal. Batteries CANNOT be replaced.

To remove the batteries before discarding product:

- Place the handle in the charging base and turn counter clockwise to remove the cap.
- Slide the battery out of the handle and remove.
- Dispose of the cell-pack in compliance with national or local regulations for environmental protection and recycling of raw materials. Do not dispose of the cell-pack by burning or incineration.

Note: During the first few days of use, your gums may bleed slightly. Any bleeding should stop within a few days. If your gums continue to bleed for a period of 2 weeks or more, discontinue use and consult your dentist.

Caution: Opening the handle will destroy the appliance and invalidate the warranty.

Note: A change in the shape of the bristles is normal after use. Toothbrush head should be replaced every 2-3 months.

ONE YEAR LIMITED WARRANTY

(Valid in USA only)

HoMedics, Inc., guarantees this product free from defects in material and workmanship for a period of one year from the date of original purchase, except as noted below.

This HoMedics product warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of HoMedics. This warranty is effective only if the product is purchased and operated in the USA. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under warranty. HoMedics shall not be responsible for any type of incidental, consequential or special damages. All implied warranties, including but not limited to those implied warranties of fitness and merchantability, are limited in the total duration of one year from the original purchase date.

To obtain warranty service on your HoMedics product, either hand deliver or mail the unit and your dated sales receipt (as proof of purchase), postpaid, along with check or money order in the amount of \$5.00 payable to HoMedics, Inc. to cover handling.

Upon receipt, HoMedics will repair or replace, as appropriate, your product and return it to you, postpaid. If it is appropriate to replace your product, HoMedics will replace the product with the same product or a comparable product at HoMedics' option. Warranty is solely through HoMedics Service Center. Service of this product by anyone other than HoMedics Service Center voids warranty.

This warranty provides you with specific legal rights. You may have additional rights which may vary from state to state. Because of individual state regulations, some of the above limitations and exclusions may not apply to you.

For more information regarding our product line in the USA, please visit: www.homedics.com

Due to continuous product improvements the pictures may differ slightly from the actual product.

© 2002 HoMedics, Inc. and its affiliated companies. All rights reserved. HoMedics® is a registered trademark of HoMedics, Inc. and its affiliated companies. PowerDent™, Quad Action™ and Smart Sensor™ are trademarks of HoMedics, Inc. and its affiliated companies. All rights reserved.

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CEPILLO DE DIENTES RECARGABLE

Mail to:
HoMedics
Service Center Dept. 168
3000 Pontiac Trail
Commerce Township, MI
48390

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DuPont® is a registered trademark of E.I. du Pont de Nemours & Company.