Automatic Wrist Blood Pressure Monitor

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TABLE OF CONTENTS

INTRODUCTORY INFORMATION:

| | |
|--|-----------|
| Important Product Notices and Safety Instructions | .3 |
| About Blood Pressure | .6 |
| What is Blood Pressure? | |
| Why Measure Your Blood Pressure? | |
| Blood Pressure Standard | .7 |
| How This Blood Pressure Monitor Works | .8 |

USAGE INFORMATION:

| | |
|--|------------|
| Name/Function of Each Part | .10 |
| Display Explanations | .11 |
| Installing Batteries | .13 |
| Applying the Cuff | .15 |
| Correct Measurement Posture | .17 |
| Measurement Procedure | .19 |
| Recalling Values from Memory | .25 |
| Clearing Values from Memory | .26 |
| Important Notes Regarding Your Blood Pressure Measurement | .27 |

| | |
|--|------------|
| Care, Maintenance & Cleaning | .29 |
| Potential for Radio/Television Interference | .31 |
| Troubleshooting | .33 |
| Specifications | .34 |
| Lifetime 5 Year Warranty | .35 |

IMPORTANT PRODUCT NOTICES AND SAFETY INSTRUCTIONS

When using your blood pressure monitor, basic precautions should always be followed. Please read and follow all instructions and warnings before using this product. Save these instructions for future reference.

- Please note that this is a home healthcare product only and it is not intended to serve as a substitute for the advice of a physician or medical professional.
- Do not use this device for diagnosis or treatment of any health problem or disease. Measurement results are for reference only. Consult a healthcare professional for interpretation of pressure measurements. Contact your physician if you have or suspect any

medical problem. Do not change your medications without the advice of your physician or healthcare professional.

- This product is not suitable for people with arrhythmias. This device may have difficulty determining the proper blood pressure for individuals with diabetes, poor circulation of blood, kidney problems or for users who have suffered from a stroke.
- Consult your physician before measuring blood pressure at the wrist if you have any of the following conditions: Severe hypertension, severe diabetes, severe arteriosclerosis, severe kidney disease or vascular disease that may compromise circulation.
- The pulse display is not suitable for checking the frequency of heart pacemakers.
- Electromagnetic interference: The device contains sensitive electronic components. Avoid strong electrical or electromagnetic

fields in the direct vicinity of the device (e.g. mobile telephones, microwave ovens). These may lead to temporary impairment of measurement accuracy.

- Use blood pressure monitor only for its intended use.
- Do not wrap the cuff around body parts other than your wrist.
- This unit is intended for use in measuring the blood pressure and pulse rate of adults. Do not use this device on infants or those who cannot express their own consent. This unit should not be operated by children.

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultatory method, within the limits prescribed by the American National Standard, *Manual, electronic, or automated sphygmomanometers*.

ABOUT BLOOD PRESSURE

What is Blood Pressure?

Blood pressure is the pressure exerted on the artery walls while blood flows through the arteries. The pressure measured when the heart contracts and sends blood out of the heart is systolic (highest) blood pressure. The pressure measured when the heart dilates with blood flowing back into the heart is called diastolic (lowest) blood pressure.

Why Measure Your Blood Pressure?

Among today's various health problems, those associated with high blood pressure are very common. High blood pressure dangerously correlates with cardiovascular diseases. Therefore, blood pressure monitoring is important for identifying those at risk.

BLOOD PRESSURE STANDARD

The US Department of Health & Human Services and the National Institutes of Health have developed a blood pressure standard that identifies areas of low- and high-risk blood pressure. This standard, however, is a general guideline as blood pressure varies among different people, age groups, etc.

Categories for Blood Pressure Levels in Adults

| | Systolic (mmHg) | Diastolic (mmHg) |
|----------------------|------------------------|-------------------------|
| Stage 2 Hypertension | ≥ 160 | ≥ 100 |
| Stage 1 Hypertension | 140~159 | 90~99 |
| Prehypertension | 120~139 | 80~89 |
| Normal | <120 | <80 |

For adults 18 and older who are not on medicine for high blood pressure, are not having a short-term serious illness, and do not have other conditions, such as diabetes and kidney disease.

To determine category of risk when systolic and diastolic readings fall into two areas, use the higher of the two numbers for classification.

There is an exception to the above definition of high blood pressure for people with diabetes and chronic kidney disease. A blood pressure of 130/80 mmHg or higher is considered high blood pressure for those individuals.

It is important that you consult with your physician regularly. Your physician will tell you your normal blood pressure range as well as the point at which you will be considered at risk. For reliable monitoring and reference of blood pressure, keeping long-term records is recommended.


HOW THIS BLOOD PRESSURE MONITOR WORKS

This monitor uses advanced technology to detect your blood pressure. With one touch of a button, the cuff will automatically inflate to block the blood flow through your artery. Next, the deflation process starts. This technology enables the monitor to automatically inflate and deflate to the appropriate level for each individual. Please note that any muscle movement during inflation or deflation will cause measurement error. When measurement is complete, the monitor will display your systolic pressure, diastolic pressure, and pulse readings.

NAME/FUNCTION OF EACH PART

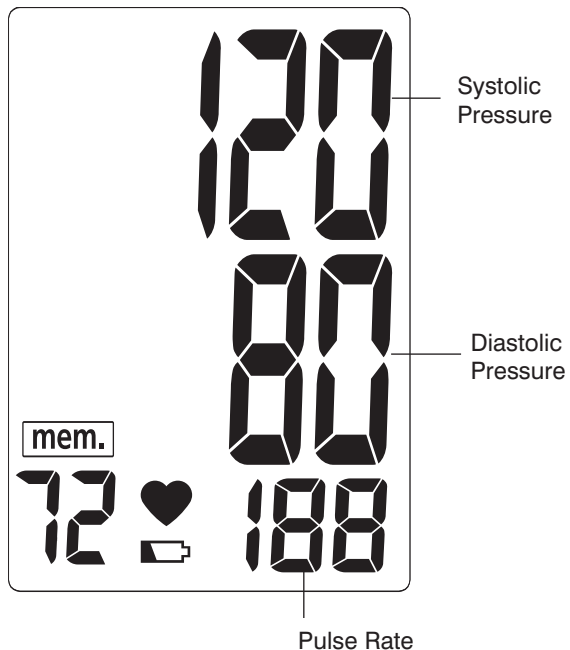DISPLAY EXPLANATIONS

Display Symbols:


| | |
|---|--|
|  | Weak Battery Symbol: Appears when batteries should be replaced. |
|  | Pulse Symbol: Shows the pulse rate per minute. |
|  | Memory Symbol: Shows the memory reading number |

If any of the following letters and numbers appear in the area that diastolic pressure should be displayed, an error has occurred with your reading. See Troubleshooting section of this manual for more information.


| | |
|---|--|
|  | Measurement Error: Measure again. Wrap the cuff correctly and keep wrist steady during measurement. |
|  | Air Circuit Abnormality: Measure again. |
|  | Pressure Exceeding 300 mmHg: Switch the unit off to clear, then measure again. |
|  | Data Error: Remove and reload the batteries. |

INSTALLING BATTERIES


1. The battery cover is located along the left side of the monitor. Place thumb and index finger on top and bottom ridges.

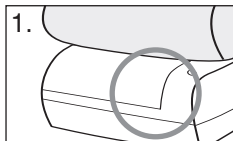2. Remove battery cover by pressing and pulling cover away from the monitor.

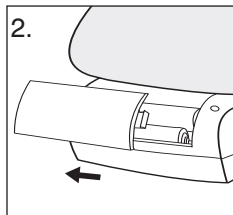
3. Install lower battery first according to the polarity indications inside the compartment. This battery should be in the bottom position as shown in figure 3.

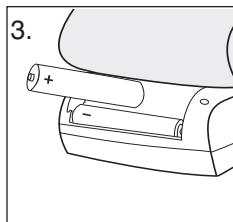
Battery Type: 2 Alkaline LR03 (AAA) size

4. Place in the other battery, according to the polarity indications inside the compartment and as shown in figure 4.

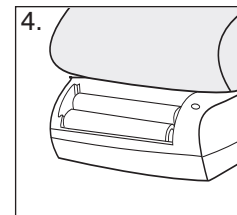5. Replace the cover by snapping it back into place.

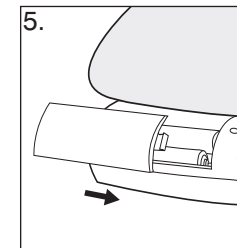Replace the batteries if:


1. The weak battery symbol appears in the display.
2. Nothing appears in the display when the power is switched on.

Replace the batteries in pairs. Remove batteries when unit is not in use for extended periods of time.


Note: Batteries are hazardous waste. Do not dispose of them together with the household garbage.

APPLYING THE CUFF

1. Remove all watches, wrist jewelry, etc. prior to attaching the wrist monitor. Clothing sleeves should be rolled up and the cuff should be wrapped around bare skin for correct measurements.
2. Apply cuff to left wrist with palm facing up as shown in Fig. A.
3. Make sure the edge of the cuff is about $1/4'' \sim 1/2''$ (1 cm) from the palm as shown in Fig. B.

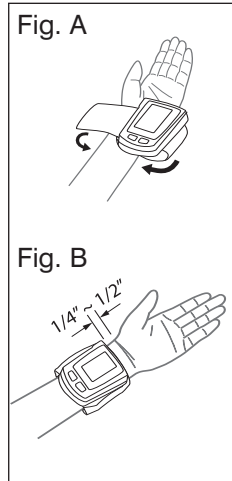4. In order to insure accurate measurements, fasten the hook and loop strap securely around your wrist so there is no extra space between the cuff and the wrist as shown in Fig. C. If the cuff is not wrapped tight enough, the measurement values will not be accurate.

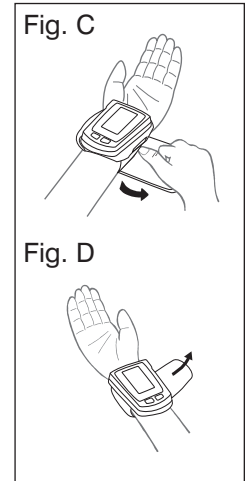5. If your physician has diagnosed you with poor circulation in your left arm, place the cuff around your right wrist as shown in Fig. D.



CORRECT MEASUREMENT POSTURE

1. Place your elbow on a table so that the cuff is at the same level as your heart as shown in Fig. A.

Note: Your heart is located slightly below your armpit.


Relax your entire body, especially between your elbow and fingers.

2. If the cuff is not at the same level as your heart or if you cannot keep your arm completely still throughout the reading, use a soft object such as a folded towel to support your arm as shown in Fig. B.

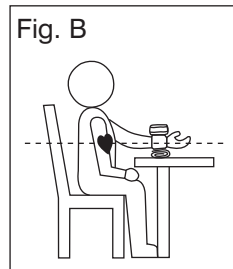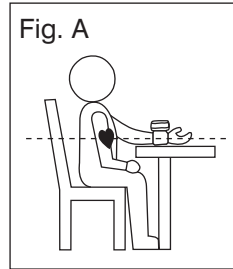3. Turn your palm upwards.
4. Sit upright in a chair, and take 5-6 deep breaths.

Avoid leaning back while the measurement is being taken as shown in Fig. C.

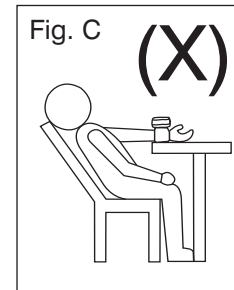
MEASUREMENT PROCEDURE

Important Notes:


- Position the monitor at the same level as your heart during measurement to ensure accurate readings.
- Blood pressure changes with every heartbeat and is in constant fluctuation throughout the day.
- Blood pressure measurement can be affected by the position of the user, his or her physiologic condition and other factors. For greatest accuracy, wait 1 hour after exercising, bathing, eating, drinking beverages with alcohol or caffeine, or smoking to measure blood pressure.

- Before measurement, it is suggested that you sit quietly for 15 minutes as measurements taken during a relaxed state will have greater accuracy. You should not be physically tired or exhausted while taking a measurement.
- During the measurement, do not talk or move your arm or hand muscles.
- See page 27 for additional notes regarding your blood pressure measurement.


1. With the cuff wrapped around your wrist, press the **START/STOP** button. Do not inflate the cuff unless it is wrapped around your wrist. All digits will light up, checking the display functions. The checking procedure will be completed after about 2 seconds.

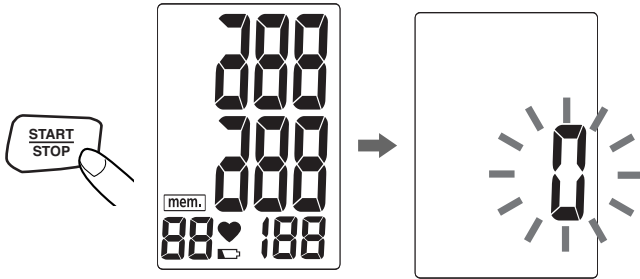
2. After all symbols disappear, the display will show a blinking "0". The monitor is "**Ready to Measure**" and will **automatically** inflate the cuff to begin measurement. The cuff will then begin deflating as measurement continues.

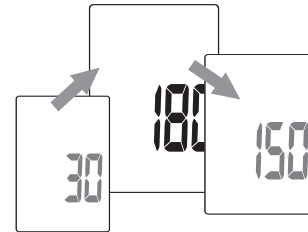Note: The monitor will re-inflate automatically if the system detects that your body requires more pressure for measurement.

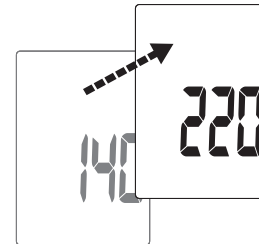3. When the measurement is completed, the cuff will deflate entirely and systolic pressure, diastolic pressure and pulse will be shown simultaneously on the LCD screen. **The measurement is then automatically stored into memory.**

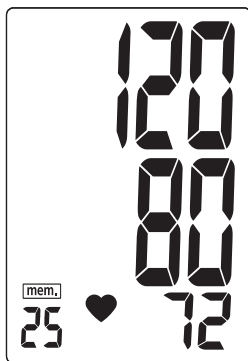Note:

1. This monitor automatically turns off approximately **1 minute** after last operation. You may also press the **START/STOP button** to turn the unit off.
2. To interrupt the measurement, you may press the **START/STOP button** (recommended) or the Memory, (**M**) button. The cuff will deflate immediately after a button is pressed.


RECALLING VALUES FROM MEMORY

This monitor can store up to 90 measurements.

1. Press the **“M” button** to access the memory.


2. Every press of the **“M” button** will recall a previous reading. The latest reading will be recalled first.

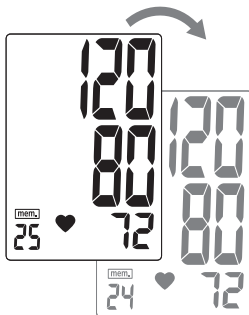CLEARING VALUES FROM MEMORY

1. Continue to press the **“M” button** for approximately 5 seconds, then the data for the selected user will be erased automatically.

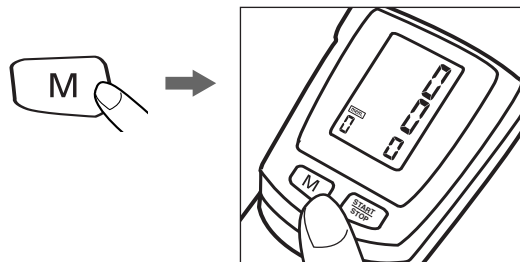Note: The data stored in memory may be erased if any of the batteries are removed. (Users should first record the data in a Blood Pressure log before replacing or removing batteries.)

IMPORTANT NOTES REGARDING YOUR BLOOD PRESSURE MEASUREMENT

- Take your reading in a comfortable environment as measurements can be affected by hot or cold temperatures. Take your blood pressure at normal body temperature.
- Do not move or talk during measurement as this can elevate readings.
- Do not move or cross legs during measurement. Keep feet flat on floor.
- Do not touch cuff or monitor during measurement procedure.
- It is suggested that you take your measurements at the same time each day and use the same arm for consistency.
- Users should wait a minimum of 5 minutes before taking additional measurements. More time may be necessary depending upon your physiology.

- The measurement results that users receive are for reference only. If users have any blood pressure concerns, please consult a physician.
- Once inflation reaches 300 mmHg, the unit will deflate automatically for safety reasons.
- This product is not suitable for people with arrhythmias.
- This device may have difficulty determining the proper blood pressure for users with diabetes, poor circulation of blood, kidney problems, or for users who have suffered a stroke.

CARE, MAINTENANCE & CLEANING

1. Clean the blood pressure monitor body and cuff carefully with a slightly damp, soft cloth. Do not press. Do not wash cuff or use chemical cleaner on it. Never use thinner, alcohol or petrol (gasoline) as cleaner.
2. Leaky batteries can damage the unit. Remove the batteries when the unit will not be used for a long time.
3. If the unit is stored near freezing, allow it to acclimate to room temperature before use.
4. This blood pressure monitor is not field serviceable. You should not use any tool to open the device nor should you attempt to adjust anything inside the device. If you have any problems with this device, please contact Heath Plus Consumer Relations (contact information on warranty page).
5. Do not immerse the unit in water as this will result in damage to the unit.

6. Do not subject the monitor and cuff to extreme temperatures, humidity, moisture, and direct sunlight. Protect from dust.
7. Do not fold the cuff tightly.
8. Do not disassemble the monitor or cuff. If in need of repair, refer to the warranty section of this manual.
9. Do not subject the monitor to extreme shocks (do not drop on floor).
10. Do not inflate the cuff unless wrapped around wrist.
11. Do not wrap the cuff around body parts other than your wrist.
12. Do not drop or insert any object into any opening.
13. This monitor may not meet its performance specifications if stored or used outside of these temperature and humidity ranges:
 - Storage Environment • Operation Environment
 - Temperature: Temperature:
14°F~140°F 50°F~104°F
(-10°C~60°C) (10°C ~40°C)
 - Humidity: Humidity:
10~90% RH max. 40~85% RH max.

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE

For USA only

This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of FCC rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation.

The product generates, uses, and can radiate radio frequency energy and, if not used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the product on and off, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and the receiver.
- Consult the dealer or an experienced radio/TV technician for help.

For Canada only

This digital apparatus does not exceed the Class B limits for radio noise emissions from digital apparatus as set out in the interference-causing equipment standard entitled “Digital Apparatus”, ICES-003 of the Canadian Department of Communications.

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TROUBLESHOOTING

If any abnormality arises during use, please check the following points:

| SYMPTOMS | POSSIBLE CAUSES | CORRECTION |
|--|--|--|
| Unit does not turn on when the START/STOP button is pushed. | Batteries have run down. | Replace them with two new AAA alkaline batteries. |
| | Battery polarities have been positioned incorrectly. | Re-insert the batteries in the correct positions. |
| EE measurement error symbol shown on display or the blood pressure value is displayed excessively low (or high). | Cuff has been placed incorrectly. | Wrap the cuff properly so that it is positioned correctly. |
| | Did you talk or move during measurement? Shaking of the wrist with the cuff on. | Measure again. Keep wrist steady during measurement. |
| E1 error symbol shown on display | Air circuit abnormality. | Measure again. |
| E2 error symbol shown on display | Inflation pressure exceeding 300 mmHg. | Switch the unit off, then measure again. |
| E3 error symbol shown on display | Data error. | Remove and reload the batteries. |

Note: If the unit still does not work, contact HealthPlus Consumer Relations. Under no circumstance should you disassemble or attempt to repair the unit by yourself. Contact information for Health Plus Consumer Relations Department can be found on the warranty page (See page 35).

SPECIFICATIONS

| | |
|---|--|
| Power Source: | Two 1.5V DC (AAA) Alkaline batteries |
| Measurement Method: | Oscillometric |
| Measurement Range: | Pressure: 40~250 millimeters Mercury (mmHg) Pulse: 40~199 beats/minute |
| Accuracy: | Pressure: ± 3 mmHg Pulse: $\pm 5\%$ of reading |
| Pressure Sensor: | Semi-conductor |
| Inflation: | Pump driven |
| Deflation: | Automatic Pressure release valve |
| Memory Capacity: | 90 memories |
| Auto-shut-off: | 1 minute after last button operation |
| Operation Environment: | Temperature: 50°F~104°F (10°C~40°C) Humidity: 40%~85% RH max. |
| Storage Environment: | Temperature: 14°F~140°F (-10°C~60°C) Humidity: 10%~90% RH max. |
| Weight: | 0.30 lbs (without batteries) |
| Wrist Circumference: | 5.3"- 8.7" (13.5 - 22cm) |
| Dimensions: | 2.6"(L) x 1.3"(W) x 3.1"(H) 67mm(L) x 32mm(W) x 78mm(H) |
|  | Type BF Device and cuff are designed to provide special protection against electrical shocks. |

Note: These specifications are subject to change without notice.

LIMITED FIVE YEAR WARRANTY

Health Plus sells its products with the intent that they are free of defects in manufacture and workmanship for a period of five years from the date of original purchase, except as noted below. Health Plus warrants that its products will be free of defects in material and workmanship under normal use and service. This warranty extends only to consumers and does not extend to Retailers.

To obtain warranty service on your Health Plus product, mail the product and your dated sales receipt (as proof of purchase), postpaid, to the following address:

Health Plus Consumer Relations
Service Center Dept. 168
3000 Pontiac Trail
Commerce Township, MI 48390

No COD's will be accepted.

Health Plus does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate Health Plus in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Health Plus.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adoption to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

For service or repair, do not return this unit to the retailer. Contact Consumer Relations at:

Mail:

**Health Plus
Service Center
Dept. 168
3000 Pontiac Trail
Commerce Township
MI 48390**

**Phone:
1-800-466-3342
Business Hours:
8:30am-5pm ET
Monday-Friday**

5 year
limited warranty

For service or repair, do not return this unit to the retailer. Contact Consumer Relations at:

Mail:

**Health Plus
Service Center
Dept. 168
3000 Pontiac Trail
Commerce
Township
MI 48390**

**Phone:
1-800-466-3342
Business Hours:
8:30am-5pm ET
Monday-Friday**

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. HEALTH PLUS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, HEALTH PLUS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS IN LIEU OF REPAIR OR REPLACEMENT.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or sales of such products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of Health Plus.

This warranty provides you with specific legal rights. You may have additional rights which may vary from country to country. Because of individual country regulations, some of the above limitations and exclusions may not apply to you.

©2007 FKA Distributing Co. and its affiliated companies, all rights reserved. Health Plus™ is a trademark of FKA Distributing Co., and its affiliated companies. All rights reserved.

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